


REVERSE BEAEBSE PEEL

PROTOCOL FOR USE

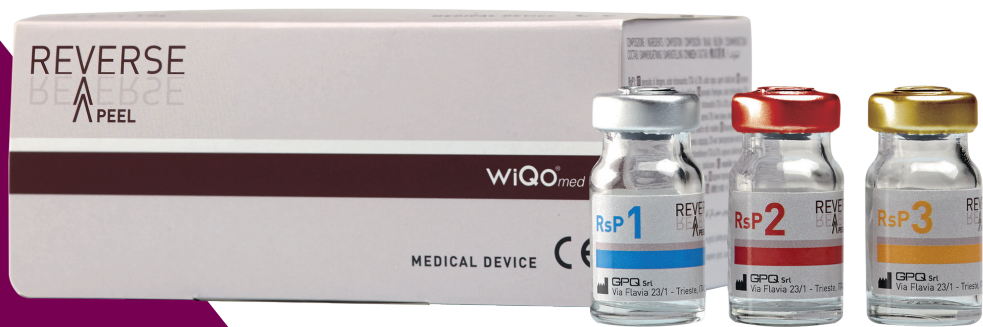

INDICATIONS

Treatment of **melasma** and **skin hyperpigmentation**.

FREQUENCY OF USE

Every 7 to 30 days, according to the doctor's assessment. Start with at least three sessions, one week apart, and then reduce the frequency, according to the intensity of melasma and its improvement.

METHOD OF USE AND APPLICATION



RsP 1

- Clean the area to be treated with WiQo **P SOLUTION**.
- Apply generously to all affected areas and massage delicately with nitrile gloves until absorbed.
- Repeat the application a second time.
- Once completely absorbed, wash the area with a soaked pad.
- Dry well.



RsP 2

- Apply superficially with the nitrile gloves on pigmented areas only, without massaging it in.
- Leave on for 3 minutes.*
- If the patient reports a burning sensation higher than 5/10 or if reddening occurs, apply the **NEUTRALISER** immediately.



N

- Shake the **NEUTRALISER** before use.
- Apply the product. Neutralisation takes place over a few seconds.
- If the patient reported burning, leave to rest until burning decreases.
- Remove with a pad soaked with water.
- Dry well.



RsP 3

- Apply superficially with the nitrile gloves on pigmented areas only, without massaging it in.
- Leave on for 3 minutes.*
- If the patient reports a burning sensation higher than 5/10 or if reddening occurs, apply the **NEUTRALISER** immediately.



N

- Shake the **NEUTRALISER** before use.
- Apply the product. Neutralisation takes place over a few seconds.
- If the patient reported burning, leave to rest until burning decreases.
- Remove with a pad soaked with water.
- Dry well.



IMMEDIATE POST-TREATMENT AND HOME TREATMENT

- Apply WiQo **REVERSE SERUM** with a delicate massage until its absorption.
- Apply WiQo **ICP CREAM** or another sun protection cream (SPF 50+).

CONTRAINDICATIONS

- Irritated skin
- Irritable skin
- Seborrhoeic dermatitis
- Atopic dermatitis
- Pregnancy and breastfeeding
- Active skin rashes and/or herpetic lesions

WARNINGS

- For medical professional use only.
- Medical device for external use.
- **The application of Phase 2 and Phase 3 must be performed with care:**
 - If redness or strong burning sensation occurs during application, immediately apply Neutraliser.
 - Do not apply Phase 3 if the skin is highly reddened or the patient reports a persistent burning sensation after the neutralisation of Phase 2.
 - If raised spots or small off-white areas appear on treated areas, neutralise immediately.
 - If the burning sensation or an intense rash persists after using the Neutraliser and subsequently rinsing the skin, apply a steroid cream.
- Do not use latex gloves as they may cause allergies. Nitrile gloves are recommended.
- Treat following a reasonable time after shaving or waxing.
- Avoid contact with eyes. In case of contact with eyes rinse thoroughly with saline solution, to be kept at hand at all time.
- Avoid contact with mucous membranes. In case of contact with mucous membranes, immediately apply the Neutraliser and then remove with a soaked pad.
- The product is corrosive.

* RsP 2 and RsP 3 can be left on for up to 15 minutes in cases of very thick skin.



**Direct line to contact dermatologists and research biologists
for scientific questions or enquiries regarding outpatient treatment:
scientificdiv.gpq@wiqo.eu**

GPQ srl

Operational Headquarters: Via Martinelli, 16 - 34015 Muggia (Trieste) - Italy
Tel.: +39 040 360408 - Fax: +39 040 761393 - gpq@wiqo.eu
www.wiqo.com